




*AeroJET-D*



*AeroJET-S*

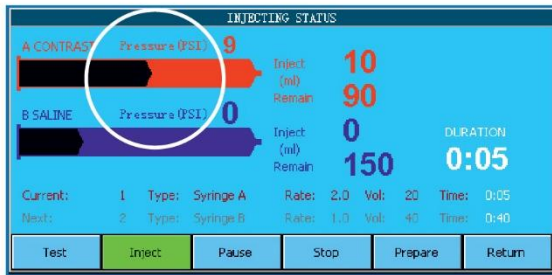


## Direct Pressure Sensor

- Real-time injection pressure display
- When the actual pressure exceeds the limit, the system will automatically stop injection and give alarm
- Effectively reduce the swelling caused by vascular rupture (See the application case)

## Local & remote real color LCD touch screen display

- Remote console: able to edit the protocols and execution of injection, synchronous display of all operations of the local console

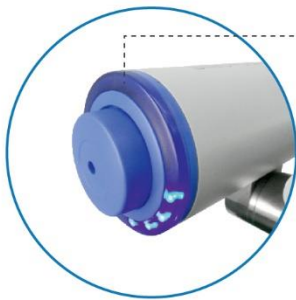

### LED rotation indicator

- LED clockwise rotation show the drug suction process
- LED counter clockwise rotation show the injection process



### Syrings heater (optional)

- Heat the contrast media and saline, making the injection process more comfortable

## Specifications of CT High-Pressure Injector

Model	AeroJET-S	AeroJET-D
Syringe	Single syringe	Dual syringe
Size of syringe	200mL sterile syringe for high-pressure injector	
Pressure limit	50 Psi~300Psi	
Flow rate	0.1mL/s~10.0mL/s Step: 0.1mL/s	
Injection volume	0.1mL~Volume of syringe	
Injection delay	0~600s	
Injection phases	1~8 phases	
Injection protocol	up to 120 protocols	
Control console	7 inch real color touch LCD for both local and remote	
Syringe heater	35°C±5°C ( ≥40min )	